


Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-09-23_15-52-30.xad

Created: 9/23/2011 3:52:30 PM
Modified: 9/23/2011 4:31:35 PM

Electrophoresis File Run Summary

[bp]

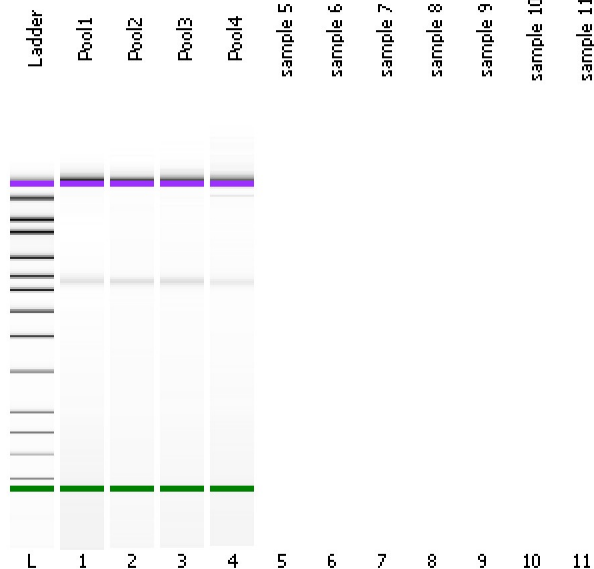

Instrument Information:





Instrument Name: Ann's baby Firmware: C.01.069
Serial#: DE24802235 Type: G2938B

Assay Information:

Assay Origin Path: D:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

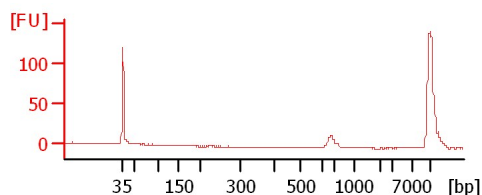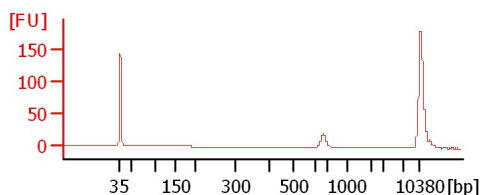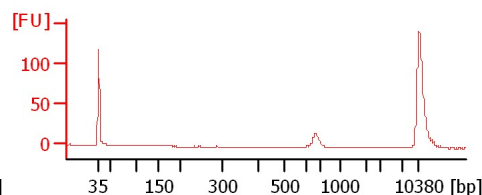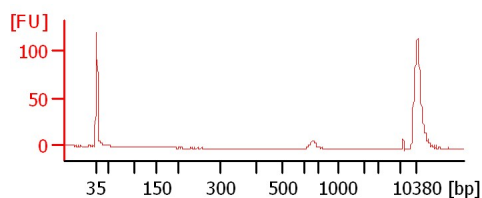Pool1**Pool2****Pool3****Pool4**

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-09-23_15-52-30.xad

Created: 9/23/2011 3:52:30 PM
Modified: 9/23/2011 4:31:35 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
Pool1		☐	✓			
Pool2		☐	✓			
Pool3		☐	✓			
Pool4		☐	✓			
sample 5		☐				
sample 6		☐				
sample 7		☐				
sample 8		☐				
sample 9		☐				
sample 10		☐				
sample 11		☐				
Ladder		☐	✓			
Chip Lot #			Reagent Kit Lot #			

Chip Comments :

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-09-23_15-52-30.xad

Created: 9/23/2011 3:52:30 PM
Modified: 9/23/2011 4:31:35 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-09-23_15-52-30.xad

Created: 9/23/2011 3:52:30 PM
 Modified: 9/23/2011 4:31:35 PM

Electropherogram Summary

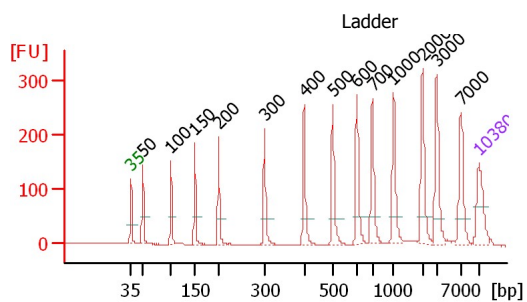

Overall Results for Ladder

Noise: 0.2

Peak table for Ladder


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-09-23_15-52-30.xad

Created: 9/23/2011 3:52:30 PM
 Modified: 9/23/2011 4:31:35 PM

Electropherogram Summary Continued ...

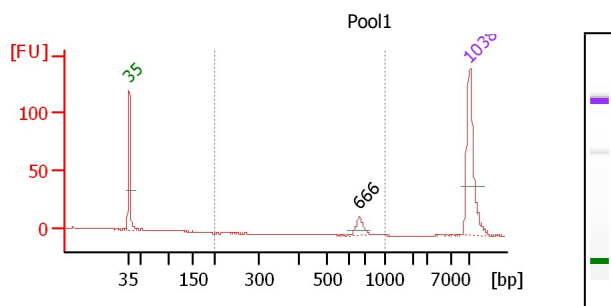

Peak table for sample 1 :

Peak	Size [bp]	Conc. [pg/μl]
1	35	125.00
2	666	15.01
3	10,380	75.00

Pool1

Molarity [pmol/l]
5,411.3

Observations

Lower Marker


Upper Marker

Overall Results for sample 1 : Pool1

Number of peaks found: 1
 Noise: 0.2
 Corr. Area 1: 25.8

Region table for sample 1 : Pool1

From To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarit y [pmol/l]	Co lo r
200 1,000	25.8	68	673	3.9	11.79	26.6	■

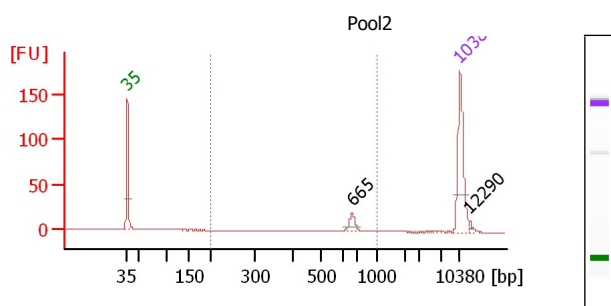

Peak table for sample 2 :

Peak	Size [bp]	Conc. [pg/μl]
1	35	125.00
2	665	16.43
3	10,380	75.00
4	12,290	0.00

Pool2

Molarity [pmol/l]
5,411.3

Observations

Lower Marker

Upper Marker

Overall Results for sample 2 : Pool2

Number of peaks found: 2
 Noise: 0.2
 Corr. Area 1: 38.3


Region table for sample 2 : Pool2

From To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarit y [pmol/l]	Co lo r
200 1,000	38.3	59	678	8.3	17.29	38.8	■

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-09-23_15-52-30.xad

Created: 9/23/2011 3:52:30 PM
 Modified: 9/23/2011 4:31:35 PM

Electropherogram Summary Continued ...

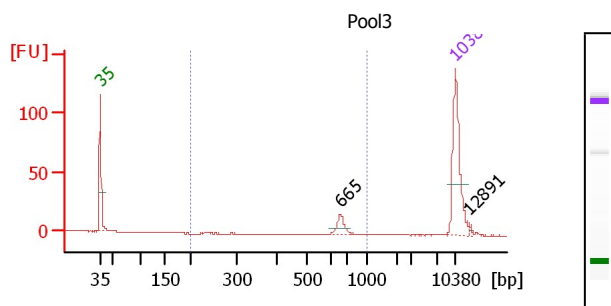

Peak table for sample 3 :


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	665	15.41	35.1	
3	10,380	75.00	10.9	Upper Marker
4	12,891	0.00	0.0	

Overall Results for sample 3 : Pool3

Number of peaks found: 2
 Noise: 0.2
 Corr. Area 1: 30.4

Region table for sample 3 : Pool3

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	30.4	69	669	5.4	15.11	34.3	Blue

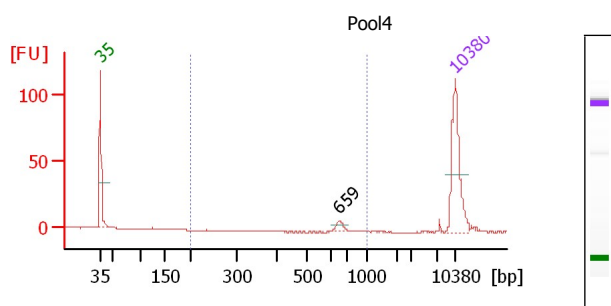

Peak table for sample 4 :

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	659	8.61	19.8	
3	10,380	75.00	10.9	Upper Marker

Overall Results for sample 4 : Pool4

Number of peaks found: 1
 Noise: 0.2
 Corr. Area 1: 13.2

Region table for sample 4 : Pool4

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	13.2	45	661	3.9	7.18	16.5	Blue

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-09-23_15-52-30.xad

Created: 9/23/2011 3:52:30 PM
Modified: 9/23/2011 4:31:35 PM

Gel Image