


Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
Modified: 10/16/2011 8:39:07 PM

Electrophoresis File Run Summary

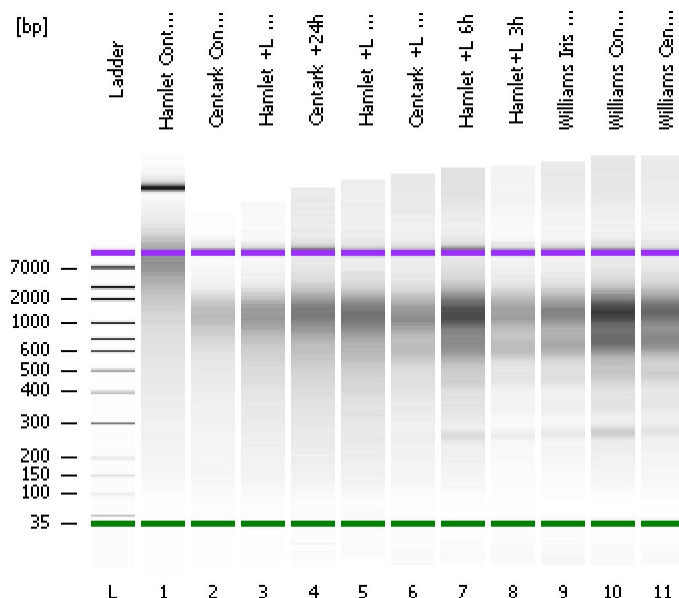

Instrument Information:


Instrument Name: Ann's baby
Serial#: DE24802235
Firmware: C.01.069
Type: G2938B

Assay Information:

Assay Origin Path: D:\Program Files\Agilent\2100 bioanalyzer\2100 expert\assays\dsDNA\High Sensitivity DNA.xsy
Assay Class: High Sensitivity DNA Assay
Version: 1.03
Assay Comments: Copyright © 2003-2010 Agilent Technologies

Chip Information:

Chip Lot #:
Reagent Kit Lot #:
Chip Comments:

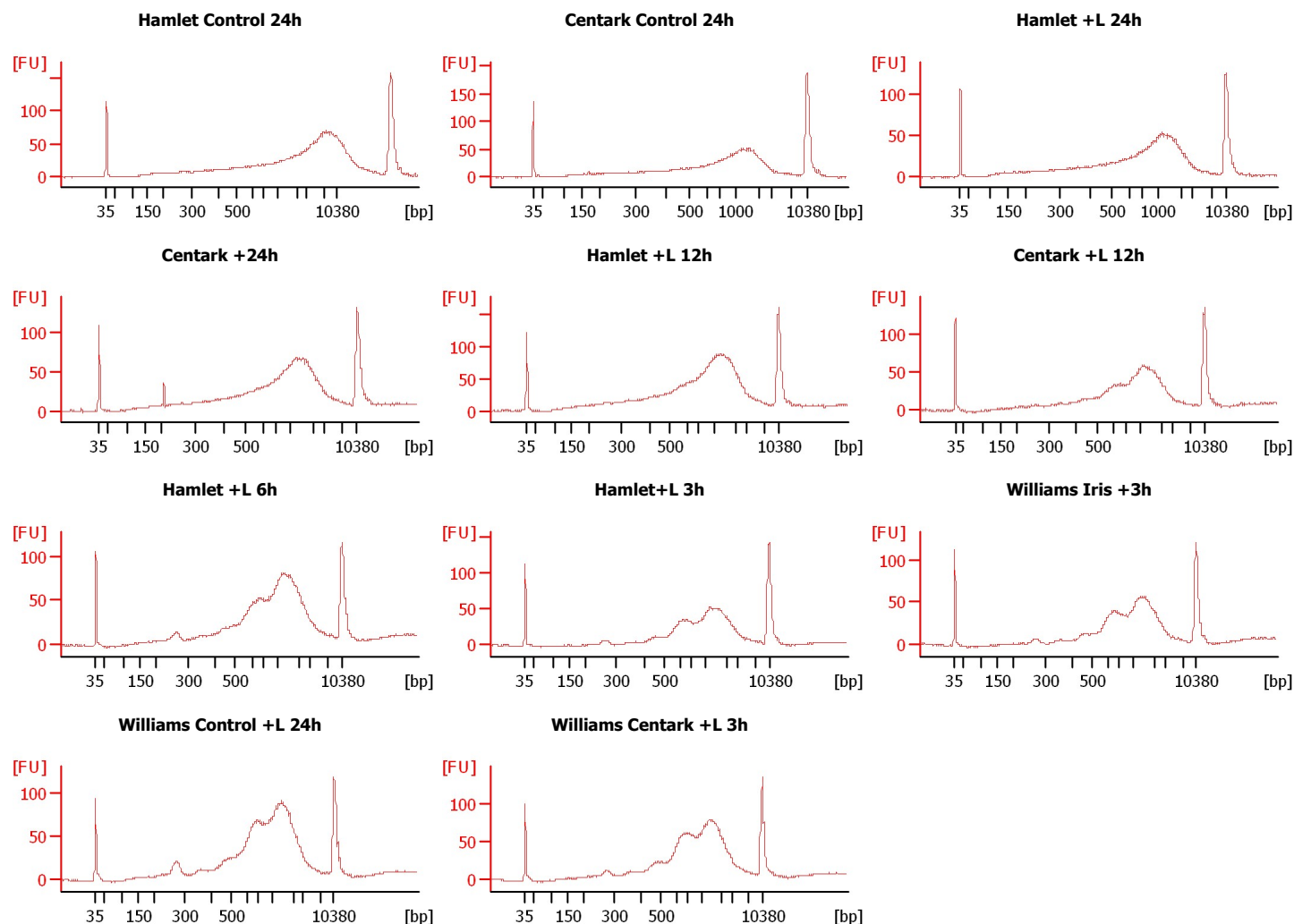

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
Modified: 10/16/2011 8:39:07 PM

Electrophoresis File Run Summary (Chip Summary)

Sample Name	Sample Comment	Rest. Digest	Status	Observation	Result Label	Result Color
Hamlet Control 24h		☐		✓		
Centark Control 24h		☐		✓		
Hamlet +L 24h		☐		✓		
Centark +24h		☐		✓		
Hamlet +L 12h		☐		✓		
Centark +L 12h		☐		✓		
Hamlet +L 6h		☐		✓		
Hamlet+L 3h		☐		✓		
Williams Iris +3h		☐		✓		
Williams Control +L 24h		☐		✓		
Williams Centark +L 3h		☐		✓		
Ladder		☐		✓		

Chip Lot #**Reagent Kit Lot #****Chip Comments :**

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
Modified: 10/16/2011 8:39:07 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/μl] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/μl] : 125
Upper Marker Concentration [pg/μl] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
 Modified: 10/16/2011 8:39:07 PM

Electropherogram Summary

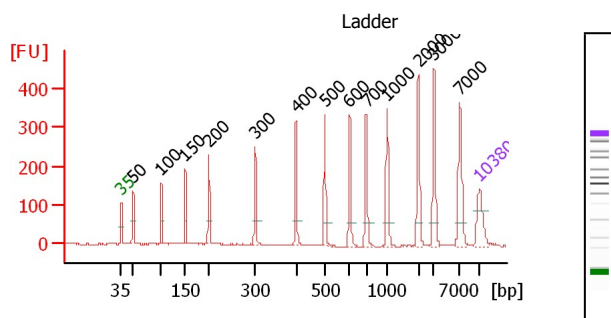


Overall Results for Ladder

Noise: 0.3

Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

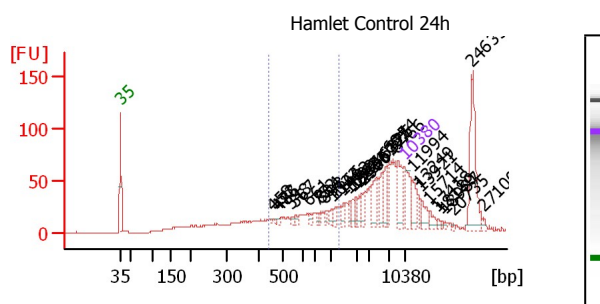

Overall Results for sample 1 : Hamlet Control 24h

Number of peaks found: 32
 Noise: 0.4
 Corr. Area 1: 300.2

Region table for sample 1 : Hamlet Control 24h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
439	1,289	300.2	19	724	31.2	378.27	886.7	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad


Created: 10/16/2011 6:53:21 PM
 Modified: 10/16/2011 8:39:07 PM

Electropherogram Summary Continued ...

Peak table for sample 1 : Hamlet Control 24h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	458	5.63	18.6	
3	473	5.69	18.2	
4	486	6.28	19.6	
5	548	9.73	26.9	
6	567	11.32	30.2	
7	661	11.84	27.1	
8	752	10.53	21.2	
9	804	10.95	20.6	
10	898	11.12	18.8	
11	1,022	17.89	26.5	
12	1,166	14.83	19.3	
13	1,462	18.42	19.1	
14	1,699	17.12	15.3	
15	1,868	18.13	14.7	
16	1,986	20.68	15.8	
17	2,224	18.51	12.6	
18	2,560	29.35	17.4	
19	4,537	70.02	23.4	
20	5,797	41.72	10.9	
21	6,325	49.13	11.8	
22	6,854	37.61	8.3	
23	8,766	80.07	13.8	
24	10,380	75.00	10.9	Upper Marker
25	11,994	0.00	0.0	

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
26	13,249	0.00	0.0	
27	13,921	0.00	0.0	
28	15,714	0.00	0.0	
29	17,418	0.00	0.0	
30	18,180	0.00	0.0	
31	18,987	0.00	0.0	
32	20,735	0.00	0.0	
33	24,635	0.00	0.0	
34	27,100	0.00	0.0	

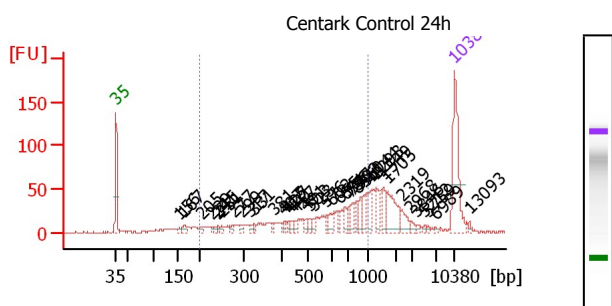

Overall Results for sample 2 : Centark Control 24h

Number of peaks found: 40
 Noise: 0.9
 Corr. Area 1: 685.2

Region table for sample 2 : Centark Control 24h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	685.2	57	575	38.2	315.95	1,081.6	Blue

Peak table for sample 2 : Centark Control 24h


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	156	5.95	57.7	
3	167	9.81	89.0	
4	205	5.08	37.6	
5	229	4.52	29.9	
6	241	4.31	27.1	
7	245	4.24	26.2	
8	251	5.34	32.2	
9	277	4.29	23.5	
10	299	5.62	28.4	
11	317	5.97	28.5	
12	331	4.41	20.2	
13	381	4.49	17.9	
14	407	4.48	16.7	
15	415	4.68	17.1	
16	427	5.17	18.3	
17	441	7.54	25.9	
18	457	5.82	19.3	
19	491	5.68	17.5	
20	503	6.78	20.4	
21	523	7.47	21.7	
22	576	6.79	17.8	
23	602	11.36	28.6	
24	655	7.22	16.7	
25	691	12.84	28.1	

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
 Modified: 10/16/2011 8:39:07 PM

Electropherogram Summary Continued ...

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
26	736	10.69	22.0	
27	777	10.79	21.0	
28	856	10.22	18.1	
29	903	11.53	19.3	
30	960	13.66	21.6	
31	1,144	18.15	24.0	
32	1,293	19.04	22.3	
33	1,549	17.19	16.8	
34	1,705	41.10	36.5	
35	2,319	13.36	8.7	
36	2,968	4.65	2.4	
37	3,846	3.30	1.3	
38	5,212	2.05	0.6	
39	5,759	1.56	0.4	
40	6,989	1.52	0.3	
41	10,380	75.00	10.9	Upper Marker
42	13,093	0.00	0.0	

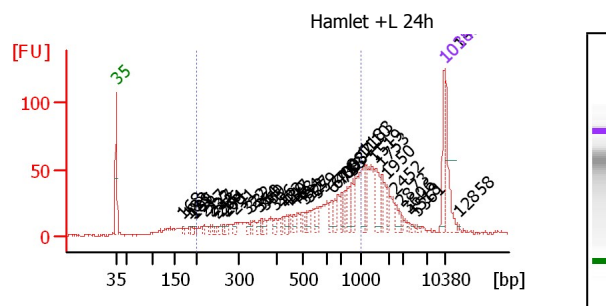

Overall Results for sample 3 : Hamlet +L 24h

Number of peaks found: 44
 Noise: 0.8
 Corr. Area 1: 673.6

Region table for sample 3 : Hamlet +L 24h

From To [bp]	Corr. Area [bp]	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	673.6	56	580	37.4	662.85	2,230.4


Peak table for sample 3 : Hamlet +L 24h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	168	6.73	60.6	
3	183	6.96	57.6	
4	189	7.19	57.5	
5	196	8.71	67.3	
6	217	7.64	53.3	
7	233	8.61	56.0	
8	237	6.92	44.3	
9	248	6.76	41.4	
10	256	8.45	50.1	
11	264	6.71	38.5	
12	271	7.17	40.1	
13	291	9.08	47.3	
14	322	11.01	51.8	
15	332	9.86	45.0	
16	354	9.81	41.9	
17	360	8.41	35.4	
18	369	8.14	33.4	
19	396	11.94	45.7	
20	422	11.46	41.2	
21	439	12.05	41.6	
22	447	13.96	47.3	
23	468	14.93	48.4	
24	494	12.82	39.3	
25	512	13.87	41.1	

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
26	527	16.95	48.7	
27	559	12.64	34.3	
28	640	23.65	56.0	
29	671	23.74	53.6	
30	709	18.52	39.5	
31	758	19.09	38.2	
32	803	29.35	55.4	
33	887	29.09	49.7	
34	1,110	32.30	44.1	
35	1,183	27.43	35.1	
36	1,519	48.35	48.2	
37	1,753	20.11	17.4	
38	1,950	17.53	13.6	
39	2,452	18.64	11.5	
40	2,873	5.02	2.6	
41	3,606	3.55	1.5	
42	4,518	2.39	0.8	
43	5,361	3.13	0.9	
44	10,148	42.48	6.3	
45	10,380	75.00	10.9	Upper Marker
46	12,858	0.00	0.0	

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
Modified: 10/16/2011 8:39:07 PM

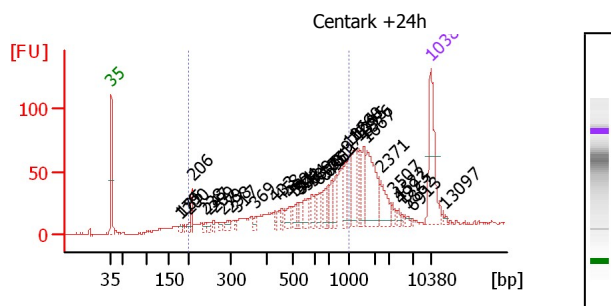Electropherogram Summary Continued ...**Overall Results for sample 4 : Centark +24h**

Number of peaks found: 37
Noise: 0.7
Corr. Area 1: 760.7


Region table for sample 4 : Centark +24h

From To [bp]	Corr. Area [bp]	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarit y [pmol/l]	Co lo r
200	1,000	760.7	56	590	36.9	467.42	1,555.1

Peak table for sample 4 :

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	179	6.02	50.9	
3	185	5.28	43.2	
4	200	7.39	56.0	
5	206	11.24	82.8	
6	238	6.05	38.4	
7	247	5.78	35.5	
8	259	5.77	33.8	
9	279	6.26	34.0	
10	298	10.87	55.2	
11	317	5.46	26.1	
12	369	7.57	31.1	
13	433	6.86	24.0	
14	458	8.05	26.6	
15	490	15.93	49.3	
16	511	12.50	37.1	
17	524	9.02	26.1	
18	542	12.75	35.7	
19	564	19.38	52.0	
20	589	12.86	33.1	
21	630	17.70	42.6	
22	670	12.56	28.4	
23	696	14.54	31.6	
24	751	20.73	41.8	
25	811	22.04	41.2	
26	975	35.88	55.8	

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
27	1,064	18.27	26.0	
28	1,260	41.32	49.7	
29	1,395	19.15	20.8	
30	1,576	29.60	28.5	
31	1,667	76.17	69.2	
32	2,371	30.75	19.7	
33	3,507	6.21	2.7	
34	4,522	3.03	1.0	
35	4,993	2.34	0.7	
36	5,862	2.63	0.7	
37	6,913	1.50	0.3	
38	10,380	75.00	10.9	Upper Marker
39	13,097	0.00	0.0	

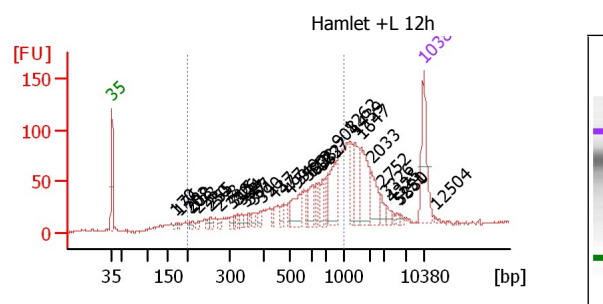**Overall Results for sample 5 : Hamlet +L 12h**

Number of peaks found: 36
Noise: 0.9
Corr. Area 1: 996.2

Region table for sample 5 : Hamlet +L 12h

From To [bp]	Corr. Area [bp]	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarit y [pmol/l]	Co lo r
200	1,000	996.2	58	597	35.9	567.45	1,832.7

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad


Created: 10/16/2011 6:53:21 PM
 Modified: 10/16/2011 8:39:07 PM

Electropherogram Summary Continued ...

Peak table for sample 5 : Hamlet +L 12h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	170	5.56	49.7	
3	176	4.90	42.1	
4	202	6.84	51.3	
5	208	6.06	44.1	
6	223	5.94	40.4	
7	250	8.81	53.3	
8	255	8.55	50.7	
9	278	10.76	58.6	
10	310	8.35	40.8	
11	315	9.05	43.5	
12	326	7.47	34.7	
13	334	9.40	42.6	
14	347	8.59	37.5	
15	357	7.43	31.6	
16	371	7.12	29.0	
17	390	12.10	47.0	
18	437	10.83	37.5	
19	477	16.73	53.1	
20	499	10.33	31.4	
21	549	47.31	130.5	
22	583	13.65	35.5	
23	612	31.20	77.3	
24	638	15.35	36.4	
25	682	22.82	50.7	
26	737	18.45	37.9	

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
27	908	73.53	122.7	
28	1,262	31.94	38.4	
29	1,439	51.27	54.0	
30	1,647	73.59	67.7	
31	2,033	39.00	29.1	
32	2,752	11.06	6.1	
33	3,326	9.24	4.2	
34	4,325	3.55	1.2	
35	5,251	3.05	0.9	
36	5,880	2.57	0.7	
37	10,380	75.00	10.9	Upper Marker
38	12,504	0.00	0.0	

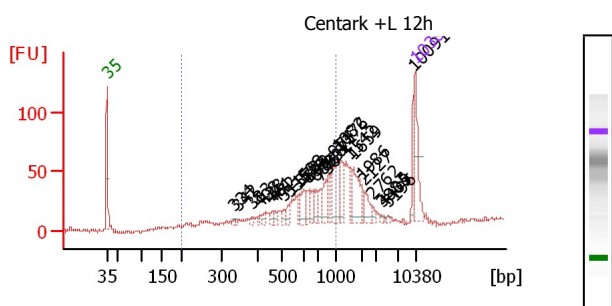

Overall Results for sample 6 : Centark +L 12h

Number of peaks found: 29
 Noise: 1.4
 Corr. Area 1: 551.2

Region table for sample 6 : Centark +L 12h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	551.2	56	630	31.8	570.16	1,669.9	Blue

Peak table for sample 6 : Centark +L 12h


Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	334	6.61	30.0	
3	342	5.42	24.0	
4	383	6.95	27.5	
5	408	6.45	24.0	
6	428	8.84	31.3	
7	461	12.11	39.8	
8	472	9.16	29.4	
9	511	10.74	31.9	
10	535	17.91	50.7	
11	578	22.52	59.0	
12	598	38.98	98.8	
13	632	25.99	62.3	
14	649	27.57	64.4	
15	689	39.74	87.4	
16	736	23.16	47.7	
17	836	41.53	75.2	
18	892	31.07	52.8	
19	974	31.24	48.6	
20	1,082	39.81	55.7	
21	1,176	46.85	60.4	
22	1,542	25.17	24.7	
23	1,659	68.15	62.2	
24	1,986	14.06	10.7	
25	2,127	28.85	20.5	

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
 Modified: 10/16/2011 8:39:07 PM

Electropherogram Summary Continued ...

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
26	2,762	5.89	3.2	
27	3,885	3.06	1.2	
28	4,409	2.15	0.7	
29	5,196	3.14	0.9	
30	10,091	38.45	5.8	
31	10,380	75.00	10.9	Upper Marker

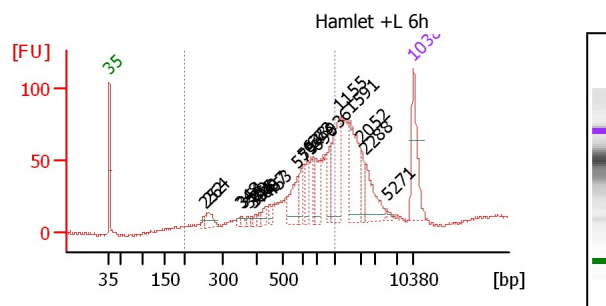

Overall Results for sample 7 : Hamlet +L 6h


Number of peaks found: 21
 Noise: 0.9
 Corr. Area 1: 730.5

Region table for sample 7 : Hamlet +L 6h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	730.5	56	636	30.5	513.77	1,476.4	Blue

Peak table for sample 7 : Hamlet +L 6h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	252	6.45	38.8	
3	264	29.10	166.9	
4	348	5.97	26.0	
5	354	7.24	31.0	
6	372	5.76	23.5	
7	384	4.54	17.9	
8	395	5.00	19.2	
9	409	6.15	22.8	
10	437	16.91	58.6	
11	453	12.21	40.8	
12	575	54.71	144.2	
13	594	23.60	60.1	
14	638	33.57	79.7	
15	673	22.33	50.3	
16	696	40.08	87.3	
17	936	37.61	60.9	
18	1,155	82.89	108.7	
19	1,591	93.67	89.2	
20	2,052	20.40	15.1	
21	2,288	56.42	37.4	
22	5,271	4.29	1.2	
23	10,380	75.00	10.9	Upper Marker

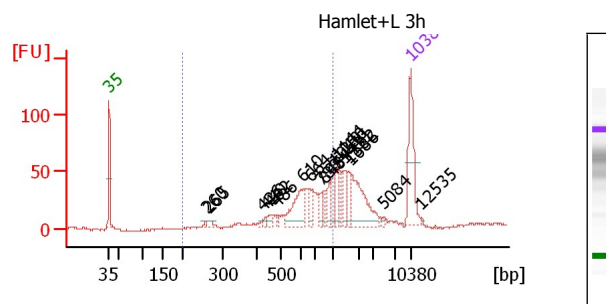

Overall Results for sample 8 : Hamlet+L 3h

Number of peaks found: 19
 Noise: 0.6
 Corr. Area 1: 449.2

Region table for sample 8 : Hamlet+L 3h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	449.2	52	662	26.2	253.24	657.4	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad


Created: 10/16/2011 6:53:21 PM
 Modified: 10/16/2011 8:39:07 PM

Electropherogram Summary Continued ...

Peak table for sample 8 :

Hamlet+L 3h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	260	5.26	30.7	
3	265	10.41	59.4	
4	426	4.02	14.3	
5	440	4.95	17.0	
6	462	9.55	31.3	
7	486	5.15	16.0	
8	610	63.56	157.8	
9	664	14.26	32.6	
10	788	15.67	30.1	
11	848	11.55	20.6	
12	954	17.08	27.1	
13	1,040	24.41	35.6	
14	1,144	17.66	23.4	
15	1,200	17.54	22.1	
16	1,431	19.98	21.2	
17	1,582	17.96	17.2	
18	1,686	80.94	72.7	
19	5,084	2.84	0.8	
20	10,380	75.00	10.9	Upper Marker
21	12,535	0.00	0.0	

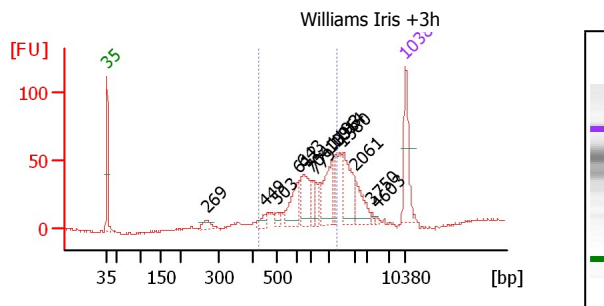

Overall Results for sample 9 : Williams Iris +3h

Number of peaks found: 14
 Noise: 0.6
 Corr. Area 1: 487.2


Region table for sample 9 : Williams Iris +3h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
429	1,326	487.2	62	796	30.1	342.26	728.2	Blue

Peak table for sample 9 :

Williams Iris +3h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	269	21.22	119.6	
3	449	15.36	51.8	
4	503	10.60	31.9	
5	612	62.97	155.8	
6	643	74.96	176.6	
7	704	23.70	51.0	
8	779	20.66	40.2	
9	1,096	68.56	94.8	
10	1,193	21.24	27.0	
11	1,314	23.26	26.8	
12	1,580	75.47	72.4	
13	2,061	43.03	31.6	
14	3,750	3.09	1.2	
15	4,603	1.86	0.6	
16	10,380	75.00	10.9	Upper Marker

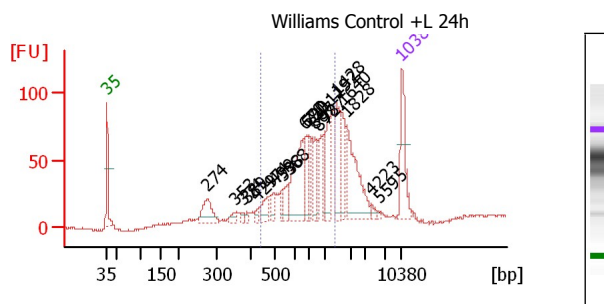

Overall Results for sample 10 : Williams Control +L 24h

Number of peaks found: 21
 Noise: 0.6
 Corr. Area 1: 893.1

Region table for sample 10 : Williams Control +L 24h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
441	1,375	893.1	57	802	30.9	617.51	1,310.7	Blue

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad


Created: 10/16/2011 6:53:21 PM
 Modified: 10/16/2011 8:39:07 PM

Electropherogram Summary Continued ...

Peak table for sample 10 :

Williams Control +L 24h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	274	69.05	382.4	
3	353	10.48	45.0	
4	374	7.30	29.5	
5	389	4.63	18.0	
6	429	7.95	28.1	
7	474	29.08	92.9	
8	499	18.02	54.8	
9	538	15.11	42.6	
10	568	31.75	84.6	
11	670	137.55	311.2	
12	690	35.25	77.4	
13	720	34.53	72.6	
14	807	38.15	71.6	
15	872	28.24	49.1	
16	987	32.67	50.1	
17	1,192	62.37	79.3	
18	1,428	71.85	76.2	
19	1,640	45.29	41.8	
20	1,828	124.42	103.1	
21	4,223	7.37	2.6	
22	5,595	2.96	0.8	
23	10,380	75.00	10.9	Upper Marker

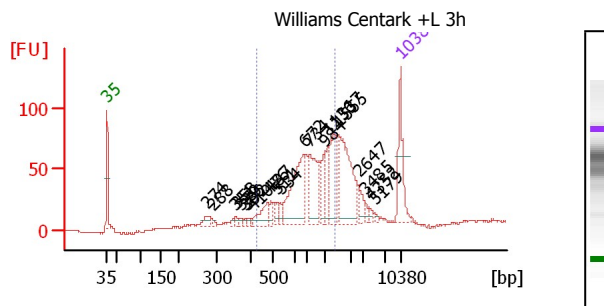

Overall Results for sample 11 : Williams Centark +L 3h

Number of peaks found: 22
 Noise: 0.8
 Corr. Area 1: 804.7

Region table for sample 11 : Williams Centark +L 3h

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
429	1,378	804.7	59	794	31.4	572.56	1,230.4	Blue

Peak table for sample 11 :

Williams Centark +L 3h

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	274	25.47	140.7	
3	288	6.61	34.8	
4	351	5.16	22.3	
5	358	8.16	34.5	
6	370	7.40	30.3	
7	386	4.86	19.1	
8	392	5.46	21.1	
9	410	4.45	16.4	
10	476	47.91	152.6	
11	487	16.53	51.4	
12	511	15.08	44.7	
13	534	11.38	32.3	
14	672	170.27	383.6	
15	734	91.67	189.2	
16	984	39.81	61.3	
17	1,153	52.11	68.5	
18	1,367	75.98	84.2	
19	1,555	157.12	153.1	
20	2,647	22.98	13.2	
21	3,485	8.10	3.5	
22	4,352	6.01	2.1	
23	5,179	3.90	1.1	
24	10,380	75.00	10.9	Upper Marker

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-10-16_18-53-22.xad

Created: 10/16/2011 6:53:21 PM
Modified: 10/16/2011 8:39:07 PM

Gel Image