

Assay Class: High Sensitivity DNA Assay
Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-01-11_15-02-07.xad

Created: 1/11/2011 3:02:07 PM
Modified: 1/11/2011 3:43:34 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 12
Minimum Visible Range [s] : 32
Maximum Visible Range [s] : 138
Start Analysis Time Range [s] : 33
End Analysis Time Range [s] : 137.5
Ladder Concentration [pg/ μ l] : 1950
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [pg/ μ l] : 125
Upper Marker Concentration [pg/ μ l] : 75
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 33.05
Integration End Time [s] : 137
Slope Threshold : 0.8
Height Threshold [FU] : 5
Area Threshold : 0.1
Width Threshold [s] : 0.6
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	35	160
2	50	210
3	100	208
4	150	221
5	200	242
6	300	270
7	400	305
8	500	306
9	600	336
10	700	321
11	1000	366
12	2000	413
13	3000	411
14	7000	400
15	10380	214

Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-01-11_15-02-07.xad

Created: 1/11/2011 3:02:07 PM
 Modified: 1/11/2011 3:43:34 PM

Electropherogram Summary

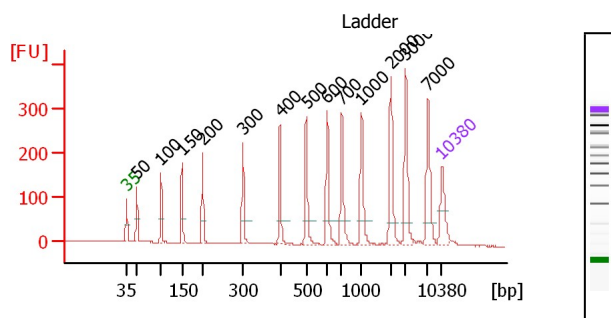


Overall Results for Ladder

Noise: 0.2

Peak table for Ladder

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	50	150.00	4,545.5	Ladder Peak
3	100	150.00	2,272.7	Ladder Peak
4	150	150.00	1,515.2	Ladder Peak
5	200	150.00	1,136.4	Ladder Peak
6	300	150.00	757.6	Ladder Peak
7	400	150.00	568.2	Ladder Peak
8	500	150.00	454.5	Ladder Peak
9	600	150.00	378.8	Ladder Peak
10	700	150.00	324.7	Ladder Peak
11	1,000	150.00	227.3	Ladder Peak
12	2,000	150.00	113.6	Ladder Peak
13	3,000	150.00	75.8	Ladder Peak

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
14	7,000	150.00	32.5	Ladder Peak
15	10,380	75.00	10.9	Upper Marker

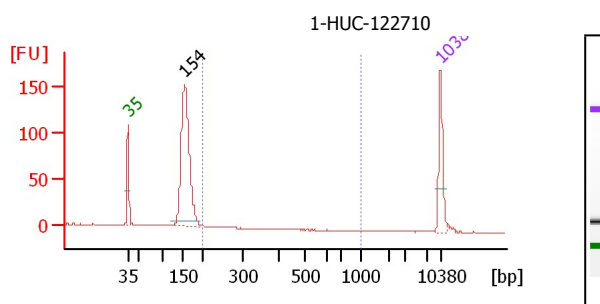

Overall Results for sample 1 : 1-HUC-122710

Number of peaks found: 1
 Noise: 0.2
 Corr. Area 1: 0.9

Region table for sample 1 : 1-HUC-122710

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	0.9	0	206	3.9	0.72	5.3	Blue


Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-01-11_15-02-07.xad

Created: 1/11/2011 3:02:07 PM
 Modified: 1/11/2011 3:43:34 PM

Electropherogram Summary Continued ...

Peak table for sample 1 : 1-HUC-122710

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.2	Lower Marker
2	154	534.28	5,246.0	
3	10,380	75.00	10.9	Upper Marker

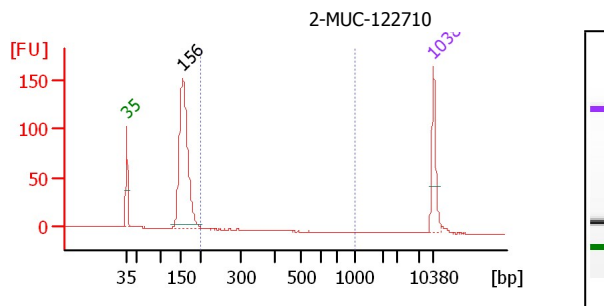

Overall Results for sample 2 : 2-MUC-122710


Number of peaks found: 1
 Noise: 0.2
 Corr. Area 1: 1.0

Region table for sample 2 : 2-MUC-122710

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	1.0	0	214	13.3	0.82	5.8	Blue

Peak table for sample 2 : 2-MUC-122710

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	156	542.37	5,267.6	
3	10,380	75.00	10.9	Upper Marker

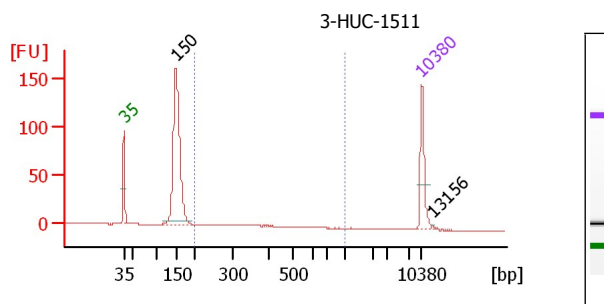

Overall Results for sample 3 : 3-HUC-1511

Number of peaks found: 2
 Noise: 0.2
 Corr. Area 1: 26.1

Region table for sample 3 : 3-HUC-1511

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	26.1	5	392	35.1	20.29	90.3	Blue


Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-01-11_15-02-07.xad

Created: 1/11/2011 3:02:07 PM
 Modified: 1/11/2011 3:43:34 PM

Electropherogram Summary Continued ...

Peak table for sample 3 : 3-HUC-1511

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	150	514.57	5,205.6	
3	10,380	75.00	10.9	Upper Marker
4	13,156	0.00	0.0	

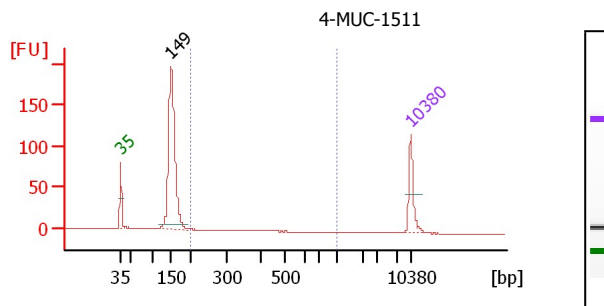

Overall Results for sample 4 : 4-MUC-1511


Number of peaks found: 1
 Noise: 0.2
 Corr. Area 1: 0.9

Region table for sample 4 : 4-MUC-1511

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	0.9	0	239	24.9	0.90	6.1	Blue

Peak table for sample 4 : 4-MUC-1511

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.2	Lower Marker
2	149	716.09	7,275.2	
3	10,380	75.00	10.9	Upper Marker

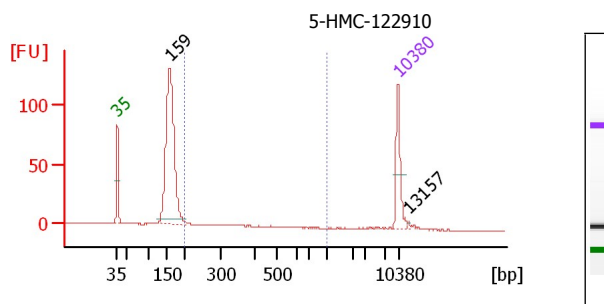

Overall Results for sample 5 : 5-HMC-122910

Number of peaks found: 2
 Noise: 0.2
 Corr. Area 1: 10.6

Region table for sample 5 : 5-HMC-122910

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	10.6	2	277	28.8	10.05	60.1	Blue


Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-01-11_15-02-07.xad

Created: 1/11/2011 3:02:07 PM
 Modified: 1/11/2011 3:43:34 PM

Electropherogram Summary Continued ...

Peak table for sample 5 : 5-HMC-122910

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	159	592.35	5,639.2	
3	10,380	75.00	10.9	Upper Marker
4	13,157	0.00	0.0	

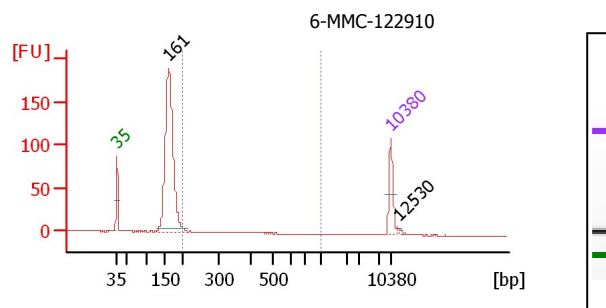

Overall Results for sample 6 : 6-MMC-122910


Number of peaks found: 2
 Noise: 0.2
 Corr. Area 1: 19.9

Region table for sample 6 : 6-MMC-122910

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	19.9	3	331	40.2	21.93	118.4	Blue

Peak table for sample 6 : 6-MMC-122910

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.2	Lower Marker
2	161	979.30	9,215.7	
3	10,380	75.00	10.9	Upper Marker
4	12,530	0.00	0.0	

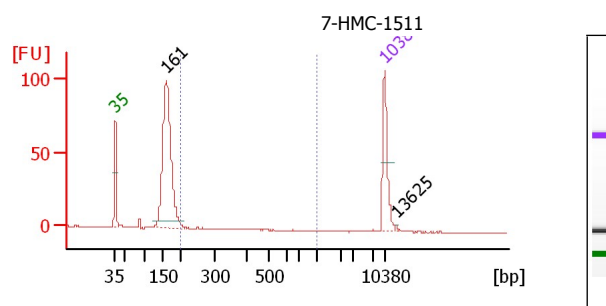

Overall Results for sample 7 : 7-HMC-1511

Number of peaks found: 2
 Noise: 0.1
 Corr. Area 1: 4.8

Region table for sample 7 : 7-HMC-1511

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	4.8	1	234	27.4	4.82	33.2	Blue


Assay Class: High Sensitivity DNA Assay
 Data Path: D:\...gh Sensitivity DNA Assay_DE24802235_2011-01-11_15-02-07.xad

Created: 1/11/2011 3:02:07 PM
 Modified: 1/11/2011 3:43:34 PM

Electropherogram Summary Continued ...

Peak table for sample 7 : 7-HMC-1511

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	161	470.17	4,427.8	
3	10,380	75.00	10.9	Upper Marker
4	13,625	0.00	0.0	

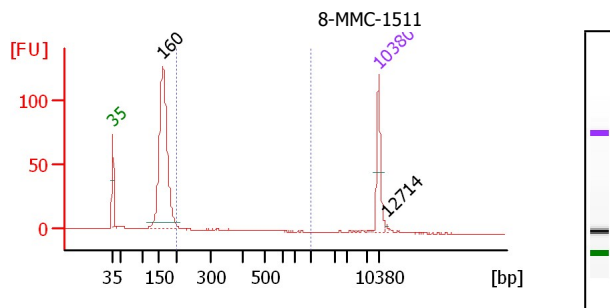

Overall Results for sample 8 : 8-MMC-1511

Number of peaks found: 2
 Noise: 0.2
 Corr. Area 1: 20.2

Region table for sample 8 : 8-MMC-1511

From [bp]	To [bp]	Corr. Area	% of Total	Average Size [bp]	Size distribution in CV [%]	Conc. [pg/μl]	Molarity [pmol/l]	Color
200	1,000	20.2	4	331	32.0	20.80	108.4	Blue

Peak table for sample 8 : 8-MMC-1511

Peak	Size [bp]	Conc. [pg/μl]	Molarity [pmol/l]	Observations
1	35	125.00	5,411.3	Lower Marker
2	160	639.49	6,058.0	
3	10,380	75.00	10.9	Upper Marker
4	12,714	0.00	0.0	