

Assay Class: DNA 7500
Data Path: D:\...-30\2100 expert_DNA 7500_DE24802235_2010-09-30_11-16-46.xad

Created: 9/30/2010 11:16:46 AM
Modified: 9/30/2010 4:10:35 PM

Electrophoresis Assay Details

General Analysis Settings

Number of Available Sample and Ladder Wells (Max.) : 13
Minimum Visible Range [s] : 20
Maximum Visible Range [s] : 94
Start Analysis Time Range [s] : 20
End Analysis Time Range [s] : 93.95
Ladder Concentration [ng/μl] : 40
Uses Standard Area for Ladder Fragments
Lower Marker Concentration [ng/μl] : 8.3
Upper Marker Concentration [ng/μl] : 4.2
Used Upper Marker for Quantitation
Standard Curve Fit is Point to Point
Show Data Aligned to Lower and Upper Marker

Integrator Settings

Integration Start Time [s] : 20
Integration End Time [s] : 93.95
Slope Threshold : 0.8
Height Threshold [FU] : 20
Area Threshold : 0.1
Width Threshold [s] : 0.5
Baseline Plateau [s] : 0.5

Filter Settings

Filter Width [s] : 0.5
Polynomial Order : 4


Ladder

Ladder Peak	Size	Area
1	50	120
2	100	47
3	300	63
4	500	81
5	700	85
6	1000	93
7	1500	101
8	2000	101
9	3000	106
10	5000	108
11	7000	109
12	10380	107

Assay Class: DNA 7500
 Data Path: D:\...-30\2100 expert_DNA 7500_DE24802235_2010-09-30_11-16-46.xad


Created: 9/30/2010 11:16:46 AM
 Modified: 9/30/2010 4:10:35 PM

Electropherogram Summary

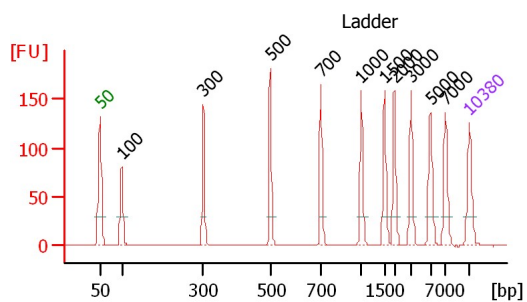

Peak table for Ladder

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	100	4.00	60.6	Ladder Peak
3	300	4.00	20.2	Ladder Peak
4	500	4.00	12.1	Ladder Peak
5	700	4.00	8.7	Ladder Peak
6	1,000	4.00	6.1	Ladder Peak
7	1,500	4.00	4.0	Ladder Peak
8	2,000	4.00	3.0	Ladder Peak
9	3,000	4.00	2.0	Ladder Peak
10	5,000	4.00	1.2	Ladder Peak
11	7,000	4.00	0.9	Ladder Peak
12	10,380	4.20	0.6	Upper Marker

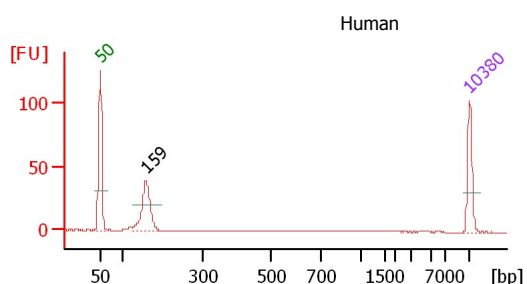


Overall Results for sample 1 : Human

Number of peaks found: 1

Peak table for sample 1 : Human

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	159	5.80	55.1	
3	10,380	4.20	0.6	Upper Marker

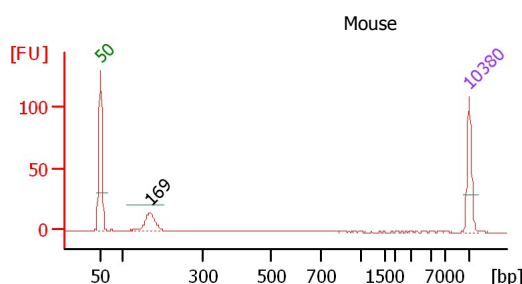

Overall Results for sample 2 : Mouse

Number of peaks found: 1

Peak table for sample 2 : Mouse

Peak	Size [bp]	Conc. [ng/μl]	Molarity [nmol/l]	Observations
1	50	8.30	251.5	Lower Marker
2	169	2.41	21.6	
3	10,380	4.20	0.6	Upper Marker

Assay Class: DNA 7500
Data Path: D:\...-30\2100 expert_DNA 7500_DE24802235_2010-09-30_11-16-46.xad

Created: 9/30/2010 11:16:46 AM
Modified: 9/30/2010 4:10:35 PM

Gel Image